

EFFICACY & SAFETY

Disposable balls for cleaning of the endoscope channels


ADVANTAGES

- High degree of cleanliness,
 Ease of use and Speed
- Reproducibility and Standardization
- Ergonomic: reduced risk of occupational diseases






medical.cbc-europe.cor

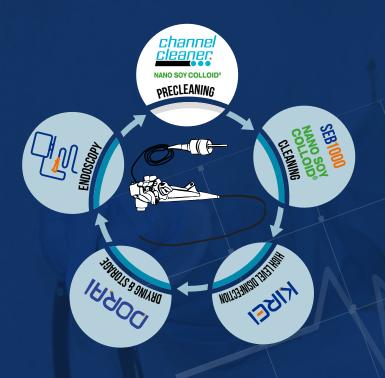

INNOVATION AND SAFETY

Manual cleaning of endoscopes is the most critical part of reprocessing. At the end of the endoscopic examination, the operator aspirates a cleaning solution through the endoscope channels to remove gross debris from patient secretions, mucus or blood.

Channel- Cleaner balls, aspirated from the distal endoscope, pass through the biopsy and aspiration channels. A trap allows the operator to verify proper ejection. The passage of the balls through the channels allows cleaning and brushing of the walls.

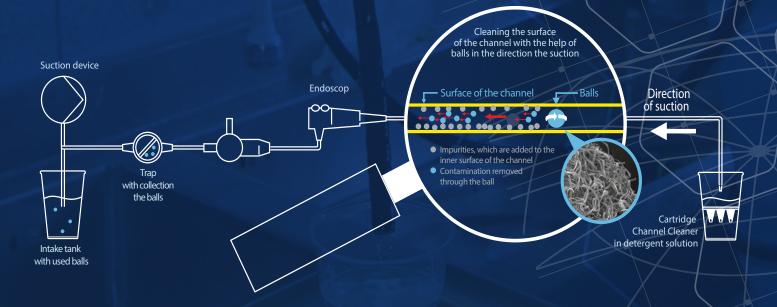

Effective, safe, repeatable and fast channel brushing, thanks to the microspheres consisting of more than a million fibers with a size of only a few microns.



CBC (Europe) GmbH

medical.cbc-europe.com

Fritz-Vomfelde-Str. 6 • 40547 Düsseldorf

info@cbc-europe.com

AVAILABLE MODELS

Model	Ball diameter
Channel Cleaner 2.0	2.0 mm
Channel Cleaner 2.8	2.8 mm
Channel Cleaner 3.2	3.2 mm

Specifications, design and accessories are subject to change without notice and without obligation on the part of the manufacturer.